

Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more tightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.

Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi



CE 2841
EN149:2001+A1:2009



1 PCS

FFP3 NR

EASY TO BREATHE



Respirator
FISH TYPE

JEDX

Functional design for uncompromised safety

Respirator FFP3 NR

With comfortable nose clip
Fish type

Instruction: 12-2020

Product number: 3805 W

JEDX

Product name: JedX 3805 W FT HLX. **Package leaflet:** Fish type 21 cm x 8,7 cm x 1,5 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtering:** solid and liquid particles; oil, dust, pollen, virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009, FFP3 NR, min 99 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi



Instruction for use:

1. Remove the respirator from the package and place it on your palm with the nose piece at the lower end of the topmost part and the elastic bands hanging down. 2. Open the respirator and place it on your face so that the nose piece fits tightly onto your nose and the lower part onto your chin. 3. Hold the respirator in place and pull both elastic bands over your head. 4. Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. 5. Shape and press the nose piece in position using your both hands to ensure a tight fit. 6. Check to make sure that the respirator fits tightly and does not leak around the edges or at the nose piece. If you detect any leaks, check the fit and re-shape the nose piece. If you are unable to make the respirator fit tightly, do not wear it.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more tightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.

Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi



CE 2841
EN149:2001+A1:2009



1 PCS

FFP2 NR

EASY TO BREATHE

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.

CE 2841
EN149:2001+A1:2009

Notified body: 2841
MNA LABORATORIES IND. TRADE.CO.LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/İstanbul | +90 216 574 07 08



JEDX

Functional design for uncompromised safety



Elastic bands
behind the head

Respirator FISH TYPE

Respirator
FFP2 NR
With comfortable nose clip
Fish type
Instruction: 12-2020
Product number: 3971 W

JEDX

Product name: JedX 3971 W FT HL. **Package leaflet:** Fish type 21 cm x 8,7 cm x 1 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtering:** solid and liquid particles: oil, dust, pollen, virus, bacterial droplets. **Filtration standard:** EN149:2001+A1:2009, FFP2 NR, min 94 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi

Instruction for use:

1. Remove the respirator from the package and place it on your palm with the nose piece at the lower end of the topmost part and the elastic bands hanging down. 2. Open the respirator and place it on your face so that the nose piece fits tightly onto your nose and the lower part onto your chin. 3. Hold the respirator in place and pull both elastic bands over your head. 4. Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. 5. Shape and press the nose piece in position using your both hands to ensure a tight fit. 6. Check to make sure that the respirator fits tightly and does not leak around the edges or at the nose piece. If you detect any leaks, check the fit and re-shape the nose piece. If you are unable to make the respirator fit tightly, do not wear it.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more rightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear a defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.



Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi

CE 2841
EN149:2001+A1:2009

Notified body: 2841

MNA LABORATORIES IND. TRADE CO. LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/İstanbul | +90 216 574 07 08

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.



CE 2841
EN149:2001+A1:2009



1 PCS

FFP2 NR

EASY TO BREATHE



Elastic bands
behind the ears

Respirator

FOLDER TYPE

JEDX

Functional design for uncompromised safety

Respirator

FFP2 NR

With comfortable nose clip

Foldable

Instruction: 12-2020

Product number: 1237 W

JEDX

Product name: JedX FFP2 NR. **Package leaflet:** Folder 15,5 x 10,5 x 0,5 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW.

Filtering: solid and liquid particles; oil, dust, pollen, virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009, FFP2 NR, min 94 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi

Instruction for use:

1. Wash your hands. Unfold the respirator and place the nose piece on your nose. 2. Pull the respirator over your cheeks and the elastic bands around your ears. 3. Adjust the nose clip to fit your nose. 4. Press the nose piece with both index fingers to prevent air leakage from the edges of the respirator.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more tightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator.
2. Do not wear defective respirator. Replace it if it is damaged in any way.
3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator.
4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections.
5. The respirator is disposable and may not be re-used.
6. When removing the mask, fold it and tie it together with the elastic bands.



Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi

CE 2841
EN149:2001+A1:2009

Notified boby: 2841

MNA LABORATORIES IND. TRADE. CO. LTD
Küçükbakkalköy Mah. Yenidoğan Cad. NO:21
Ataşehir/İstanbul | +90 216 574 07 08

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.



CE 2841
EN149:2001+A1:2009



EASY TO BREATHE

1 PCS

FFP3 NR

Respirator

FOLDER TYPE



Respirator

FFP3 NR

With comfortable nose clip

Foldable

Instruction: 12-2020

Product number: 1523 W



Instruction for use:

Product name: JedX FFP3 NR. **Package leaflet:** Folder 15,5 x 10,5 x 1,5 cm with head clip, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtering:** solid and liquid particles: oil, dust, pollen, virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009, FFP3 NR, min 99 %. The declaration of conformity in accordance with the Personal Protection, Regulation can be found at www.jedx.fi

1. Wash your hands. **2.** Unfold the respirator, place the nose piece on your nose. **3.** Pull the respirator over your cheeks and pull the elastic bands around your ears. **4.** Adjust the nose clip to fit your nose. **5.** Press the nose piece with both index fingers to prevent air leakage from the edges of the respirator.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more rightly on your face, if you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. **2.** Do not wear defective respirator. Replace it if it is damaged in any way. **3.** When wearing the respirator, stay away from open fire to avoid deformation of the respirator. **4.** If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. **5.** The respirator is disposable and may not be re-used. **6.** When removing the mask, fold it and tie it together with the elastic bands.

Warning The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi




Notified boby: 2841

MNA LABORATORIES IND. TRADE. CO. LTD
Küçükbakkalköy Mah. Yenidoğan Cad. NO:21
Ataşehir/İstanbul | +90 216 574 07 08

CE 2841
EN149:2001+A1:2009

1 PCS

FFP3 NR

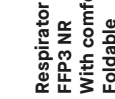

Instruction: 12-2020
Product number: 1609 W

Instruction for use: **1.** Wash your hands. Unfold the respirator, place the nose piece on your nose. **2.** Pull the respirator over your cheeks and pull both elastic bands around your head. **3.** Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. **4.** Adjust the nose piece to fit your nose. **5.** Press the nose piece with both index fingers to prevent air leakage from the edges of the respirator.



Product name: JedX 1609 W HL.V. **Package leaflet:** Folder 15,5 x 10,5 x 1,5 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtration:** solid and liquid particles; oil, dust, pollen, virus, bacterial droplets. **Filtration standard:** EN149:2001+A1:2009, FFPS NR, min 99 %. The declaration of conformity in accordance with the Personal Protection Regulation can be found at www.jedx.fi

Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more rightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.

Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi



Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.

Notified body: 2841

MNA LABORATORIES IND. TRADE. CO. LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/İstanbul | +90 216 574 07 08

CE 2841
EN149:2001+A1:2009



CE 2841
EN149:2001+A1:2009



1 PCS

FFP3 NR

EASY TO BREATHE



Elastic bands
behind the head

Respirator FOLDER TYPE

JEDX

Functional design for uncompromised safety

Respirator
FFP3 NR
With comfortable nose clip
Foldable

Instruction: 12-2020
Product number: 1625 W

JEDX

Product name: JedX 1625 W HL. **Package leaflet:** Folder 15,5 x 10,5 x 0,5 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtering:** solid and liquid particles; oil, dust, pollen, virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009, FFP3 NR, min 99 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi



Instruction for use: 1. Wash your hands. Unfold the respirator, place the nose piece on your nose. 2. Pull the respirator over your cheeks and pull both elastic bands around your head. 3. Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. 4. Adjust the nose piece to fit your nose. 5. Press the nose piece with both index fingers to prevent air leakage from the edges of the respirator.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more rightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.



Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi

CE 2841
EN149:2001+A1:2009

Notified body: 2841

MNA LABORATORIES IND. TRADE. CO. LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/İstanbul | +90 216 574 07 08

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.



CE 2841
EN149:2001+A1:2009



1 PCS

FFP2 NR

EASY TO BREATHE



Respirator FOLDER TYPE

JEDX

Functional design for uncompromised safety

Respirator
FFP2 NR
With comfortable nose clip
Foldable

Instruction: 12-2020
Product number: 1758 W

JEDX

Product name: JedX 1758 W HLV. **Package leaflet:** Folder 15,5 x 10,5 x 1,5 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtering:** solid and liquid particles; oil, dust, pollen, virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009, FFP2 NR, min 94 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi



Instruction for use: 1. Wash your hands. Unfold the respirator, place the nose piece on your nose. 2. Pull the respirator over your cheeks and pull both elastic bands around your head. 3. Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. 4. Adjust the nose piece to fit your nose. 5. Press the nose piece with both index fingers to prevent air leakage from the edges of the respirator.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more rightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.

Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi



EASY TO BREATHE

1 PCS

FFP2 NR

CE 2841
EN149:2001+A1:2009



Respirator
FFP2 NR
With comfortable nose clip
Foldable

Instruction: 12-2020
Product number: 1852 W

JEDX

Product name: JedX 1852 W HL. **Package leaflet:** Folder 15,5 x 10,5 x 0,5 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtering:** solid and liquid particles; oil, dust, pollen, virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009, FFP2 NR, min 94 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.

Notified body: 2841

MNA LABORATORIES IND. TRADE. CO. LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/Istanbul | +90 216 574 07 08

CE 2841
EN149:2001+A1:2009



Elastic bands
behind the head

JEDX
Functional design for uncompromised safety



Instruction for use: 1. Wash your hands. Unfold the respirator, place the nose piece on your nose. 2. Pull the respirator over your cheeks and pull both elastic bands around your head. 3. Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. 4. Adjust the nose piece to fit your nose. 5. Press the nose piece with both index fingers to prevent air leakage from the edges of the respirator.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more rightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear a defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.



Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi

CE 2841
EN149:2001+A1:2009

Notified body: 2841

MNA LABORATORIES IND. TRADE. CO. LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/İstanbul | +90 216 574 07 08

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.



CE 2841
EN149:2001+A1:2009



1 PCS

FFP2 NR

EASY TO BREATHE

Respirator FOLDER TYPE



JEDX

Functional design for uncompromised safety

Respirator

FFP2 NR

With comfortable nose clip

Foldable

Instruction: 12-2020

Product number: 2793 W

JEDX

Product name: JedX FFP2 NR SPORT W. **Package leaflet:** Folder
15,5 x 10,5 x 0,5 cm, 1 pcs / bag. **Material:** Spunbond NW,
Meltblown NW. **Filtering:** solid and liquid particles; oil, dust, pollen,
virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009,
FFP2 NR, min 94 %. The declaration of conformity in accordance
with the Personal Protection. Regulation can be found at www.jedx.fi

Instruction for use:

1. Wash your hands. 2. Unfold the respirator and place the nose piece on your nose. 3. Pull the respirator over your cheeks and the elastic bands around your ears. 4. Adjust the nose clip to fit your nose. 5. Press the nose piece with both index fingers to prevent air leakage from the edges of the respirator.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more tightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.

Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi



CE 2841
EN149:2001+A1:2009

Notified body: 2841
MNA LABORATORIES IND. TRADE.CO.LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/İstanbul | +90 216 574 07 08



CE 2841
EN149:2001+A1:2009



1 PCS

FFP2 NR

EASY TO BREATHE



Respirator
FISH TYPE

JEDX

Functional design for uncompromised safety

Respirator

FFP2 NR

With comfortable nose clip


Fish type

Instruction: 12-2020

Product number: 3614 W

JEDX

Product name: JedX 3614 W FT HLX. **Package leaflet:** Fish type 21 cm x 8,7 cm x 1,5 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtering:** solid and liquid particles; oil, dust, pollen, virus, bacterial, droplets. **Filtration standard:** EN149:2001+A1:2009, FFP2 NR, min 94 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi

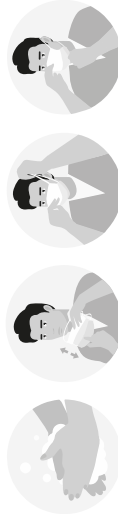

Instruction for use:

1. Remove the respirator from the package and place it on your palm with the nose piece at the lower end of the topmost part and the elastic bands hanging down. 2. Open the respirator and place it on your face so that the nose piece fits tightly onto your nose and the lower part onto your chin. 3. Hold the respirator in place and pull both elastic bands over your head. 4. Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. 5. Shape and press the nose piece into position using your both hands to ensure a tight fit. 6. Check to make sure that the respirator fits tightly and does not leak around the edges or at the nose piece. If you detect any leaks, check the fit and re-shape the nose piece. If you are unable to make the respirator fit tightly, do not wear it.



Attention

1. Place the respirator on your face as instructed and take a deep breath to check that the respirator does not leak around the edges. If it leaks, place it more tightly on your face. If you are unable to fit the respirator tightly on your face, you will not be able to protect yourself from airborne impurities with this respirator. 2. Do not wear defective respirator. Replace it if it is damaged in any way. 3. When wearing the respirator, stay away from open fire to avoid deformation of the respirator. 4. If you must remove the respirator temporarily, avoid touching its outer surface to avoid infections. 5. The respirator is disposable and may not be re-used. 6. When removing the mask, fold it and tie it together with the elastic bands.

Manufacturer:

JedX Medcare, Köysikuja 1,
01640 Vantaa, Finland
+358 20 741 5180 | www.jedx.fi



EASY TO BREATHE

1 PCS

FFP3 NR

CE 2841
EN149:2001+A1:2009



Respirator FFP3 NR

With comfortable nose clip

Fish type

Instruction: 12-2020

Product number: 3972 W

JEDX

Product name: JedX 3972 W FT HL. **Package leaflet:** Fish type 21 cm x 8,7 cm x 1 cm, 1 pcs / bag. **Material:** Spunbond NW, Meltblown NW. **Filtration:** solid and liquid particles: oil, dust, pollen, virus, bacterial droplets. **Filtration standard:** EN149:2001+A1:2009, FFP3 NR, min 99 %. The declaration of conformity in accordance with the Personal Protection. Regulation can be found at www.jedx.fi

Warning

The respirator only filters the impurities listed in the instructions. The product provides no protection against gases, vapors, or oxygen shortage. Misuse may result in illness or even death. Maximum period of use 8 hours. Shelf life 5 years for date of manufacture when stored correctly. Store in a dry (RH 0...80%, -30...+40°C) well-ventilated location protected from sunlight.

Use-by date: See package.

CE 2841
EN149:2001+A1:2009

Notified body: 2841
MNA LABORATORIES IND. TRADE.CO.LTD
Küçükbakkalköy Mah. Yenidoğan Cad.NO:21
Ataşehir/İstanbul | +90 216 574 07 08

JEDX

Functional design for uncompromised safety



Respirator FISH TYPE



Instruction for use:

1. Remove the respirator from the package and place it on your palm with the nose piece at the lower end of the topmost part and the elastic bands hanging down. 2. Open the respirator and place it on your face so that the nose piece fits tightly onto your nose and the lower part onto your chin. 3. Hold the respirator in place and pull both elastic bands over your head. 4. Place the elastic bands so that the upper band sits on your head, above your ears and the lower band is on the back of your head below your ears. 5. Shape and press the nose piece into position using your both hands to ensure a tight fit. 6. Check to make sure that the respirator fits tightly and does not leak around the edges or at the nose piece. If you detect any leaks, check the fit and re-shape the nose piece. If you are unable to make the respirator fit tightly, do not wear it.

